

Product Brief

ConnectedHealth™ Health Manager Application

ConnectedHealth's *Health Manager Application* ("HMA") automates the collection of medical information by patients at home.



The HMA captures data from medical devices such as Glucose Meters, Weighing Scales, Blood Pressure Monitors, and Pulse Oximeters, via Bluetooth™ technology. The HMA supports both Continua™ compliant devices as well as proprietary device communication protocols

HMA runs on Android-powered mobile phones and tablets that support Bluetooth, and Bluetooth Health Device Profile (for Continua™ devices).

The HMA runs as an Android background service allowing, not requiring, any user intervention to capture data. The user simply stands on a weighing scale, takes a blood pressure or Glucose reading and the data is seamlessly transferred to the ConnectedHealth™ eHealth Platform, or any Continua™ compliant WAN Receiver framework.

Specific versions of the HMA also support communication with Qualcomm Life's 2Net Platform. The HMA incorporates security features; supporting the encryption of on-device data and encryption of communication path.

Using the Health Manager Application:

- A patient can automatically capture data from medical devices, such as Glucose Meters, Weighing Scales, Blood Pressure Monitors, etc. within a home setting; and transfer the data to the ConnectedHealth™ eHealth Platform or Qualcomm Life's 2Net Platform.
- Clinicians can extract and review historical data utilizing a standalone ConnectedHealth™ Dashboard application, or directly through secure webpages on ConnectedHealth™ eHealth Platform. Patients can review their own personal data through secure web pages on the same platform.



Health Manager Application Brief

- Automates collection of medical data by patients at home using Continua™ and medical devices using proprietary communication protocols.
- Automatically uploads data to ConnectedHealth eHealth Platform.
- Data is encrypted on device and during transmission.
- Supports Bluetooth™ Serial Port Profile and Health Device Profiles spec.
- Supports Continua™ IEEE 11073 protocol, and Continua WAN sender.

SINGAPORE ConnectedHealth


Unit 03-20/22
67 Ayer Rajah Crescent
Singapore 139950
Tel: +65 6777 9750
Fax: +65 6778 9005
Email: sales@connhealth.com
Web: www.connhealth.com

USA Sovran Inc.

3664 Grand Ave, Suite A
Oakland, CA 94610
USA
Tel: +1 510 645 9377
Fax: +1 510 645 9379
Email: sales@sovrans.net

EUROPE ARS Software GmbH

Starnberger Str. 22
D-82131 Gauting
Munich, Germany
Tel: +49 (89) 893 4130
Fax: +49 (89) 8934 1310
Email: info@ars2000.com

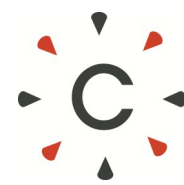


Product Brief

Supported Wireless Medical Devices

Blood Pressure Monitors	A&D UA-767 PBT-C** A&D UA-767 PBT Omron BP792IT (US), 708-BT (EU) FORA D40 2-in-1 BPM and Glucose meter**
Weighing Scales	A&D Medical UC-321PBT-C** A&D Medical UC-321PBT Omron HBF-206IT (US), BF-206BT (EU)
Pulse Oximeters	Nonin Onyx II 9560**
Glucose Meters	Entra MyGlucoHealth

Please contact ConnectedHealth™ for the latest list of supported devices ** = Continua™ compliant devices

How it all comes together

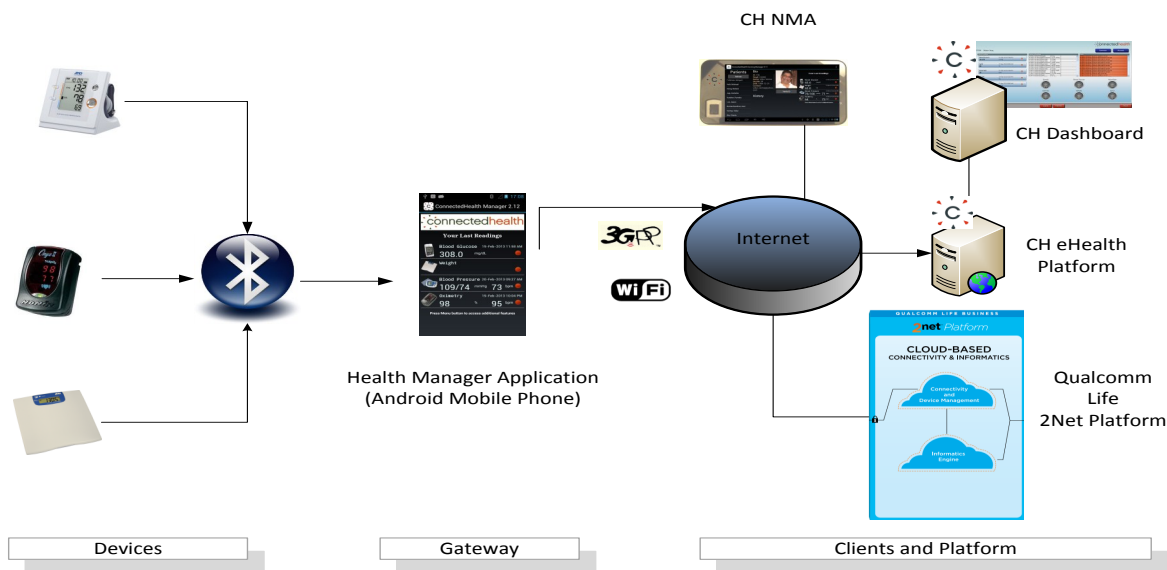

Part Numbers

Part Number	Description
CH-411-02	Android Health Manager (Proprietary Protocols)
CH-421-02	Android Health Manager (Continua™ Compliant)
CH-401-01	Android 2 Net HMA

All trademarks and trade names are those of their respective owners. ConnectedHealth makes no warranties with respect to this document and disclaims any implied warranties of merchantability and fitness for a particular purpose. Information in this document is subject to change without notice. While every attempt is made to ensure accuracy, ConnectedHealth shall not be liable for any errors or omissions that may appear in this document.